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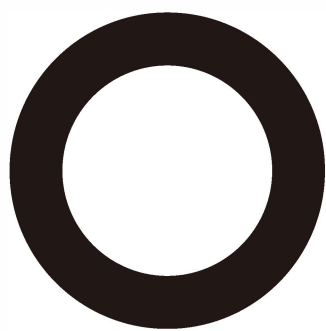
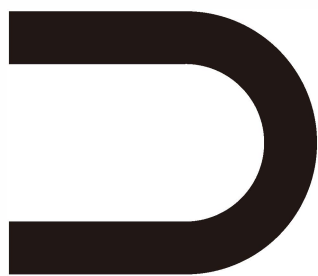
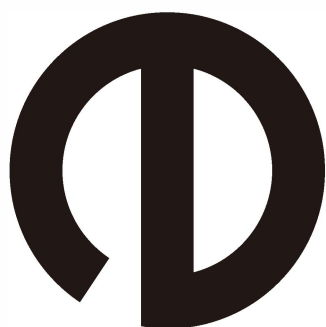
Производитель: Shenzhen Yunding Information Technology Co., Ltd.
Manufacturer Address: 28G, Building 3, Dachong Business Center (Phase 3), No.18, Dachong 1st Road, Dachong Community, Yuehai Street, Nanshan District, Shenzhen City, Guangdong Province, China

ENGLISH

Thank you for using the Oclean electric toothbrush

Important safety tips

1. Children over 8 years of age and above and persons with physical, sensory or intellectual defects, or without experience and knowledge may use the appliance under the supervision or guidance of the safe use and aware of the hazards. Children are not allowed to play with this appliance. Children shall not perform cleaning and user maintenance without supervision.
2. This product can only use Oclean special brush head, will not bear all the consequences of non-certified third party brush head.
3. This product is only personal cleaning products and is not recommended for medical use.
4. This product is only for indoor use, do not be placed in the direct sunlight area, it is recommended to be put in the room temperature dry ventilation place.
5. Please keep the body dry before charging, dry with soft cloth, and clean the charging interface with dry cloth regularly to prevent stains from deposition.
6. The appliance can only be supplied at the safety ultra-low voltage corresponding to the sign on the appliance.
7. Do not expose the device or battery to too high temperature. The normal temperature range of battery charging is 0 to 45 degrees.
8. Pay attention to the battery power supply equipment or battery terminals do not be touched by metal objects to avoid the risk of short circuit

Brush head

When replacing the brush head, please pull out or insert vertically, do not rotate left or right

Dentists recommend changing the brush head once every 2 to 3 months

Oral sensitive patients, such as the following circumstances should pay attention to

People who have just had oral and periodontal surgery or are undergoing treatment.

People with other symptoms in the mouth, including weak dentures or persistent massive bleeding from their gums.

Other oral cavity sensitive and vulnerable populations, etc.

Use of this product is not recommended, or to consult a professional dental practitioner before use.

This product is not suitable for the disabled, mentally ill or lack of relevant experience [including children (0-12 years old)].

Unless used under the supervision of their guardian.

Use Guidance

1. Click the button to turn on and off the electric toothbrush. Click the key 5S again to select the brushing mode, the screen and the corresponding mode symbols will follow the switch. In the brushing state, long press button 2S to switch mode.

2.5 brushing modes:

"sensitive care" -Designed for people with sensitive teeth, gentle cleaning while reducing irritation to sensitive areas.

"Quick Cleaning" -Suitable for daily efficient cleaning, effective removal of food residue and dirt.

"Deep cleaning" -Go deep between the teeth and hard-to-reach areas to forcefully remove stubborn stains and bacteria.

"Bright white polishing" -Strong cleaning makes the tooth surface more smooth and bright white.

"Private customization" -According to their own gums and oral status, users can intelligently develop oral care plans on the APP side, such as: brushing intensity and brushing duration setting, etc. After completing the setting, the toothbrush side will be synchronized with the current set mode.

Therefore, you can choose different models according to your preference and the advice of the dentist.

3. Wifi light breathe when the APP conducts the Wifi network. Wifi in the connection, the white light breathing, Wifi connection success, the white light is always on.

4. Restore the factory Settings: long press the button for 8 seconds to restore the factory Settings of the toothbrush.

Google Play Store & iOS APP Store

Install the APP: Scan the QR code to download the "Oclean Care+" APP, please search for and connect the toothbrush in the APP; Note: The connection is invalid in the mobile phone -system Settings-Bluetooth Settings:-

APP Connection (Home page-Add Device-Select Wifi-Enter password-Connect in-Connect successfully)

Brushing your teeth healthily

- Move the bristles at a 45 degree Angle and slowly move the toothbrush;
- During the process of brushing, the vibration effect of the toothbrush will have a short pause every 30 seconds, indicating that the brushing position should be changed;
- If there is no active pause, it will automatically turn down after brushing mode.

Charging method:

Before charging, wipe the bottom of the toothbrush and open the Type C port of the bottom of the machine to charge. When charging, the screen interface shows the charging state at the same time, and it can be filled within 4 hours.

Low power prompt:

When the toothbrush is less than 20%, the screen reminds you to charge.

1. Do NOT dismantle, impact, squeeze or put into the fire.
2. If the product is seriously bulging, do not continue to use it.
3. Do NOT expose electrical appliances or batteries to too high temperatures.
4. The battery is prohibited after immersion.
5. Pay attention to the danger of short circuit of the terminal of the battery power supply appliance or the battery by the metal object.

FCC Warning

This device complies with part I of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part I of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment *off* and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

CE Warning

1.the appliance is not to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction

2.children being supervised not to play with the appliance

3.Signal words WARNING, CAUTION, DANGER in uppercase having a height as specified : Uppercase letter of the text explaining the signal word not smaller than 1,6 mm : Moulded in, engraved, or stamped markings either raised above or have a depth below the surface of at least 0,25 mm, unless

4.This appliance contains batteries that are non-replaceable

Before the product is discarded, the battery must be removed from the product in the following steps: 1. Open the bottom cover waterproof cover, remove the screw plug.

2. Unscrew the bottom screw with a screwdriver and open the bottom cover.

3. Use a screwdriver to open the buckle on both sides, remove the movement, cut the battery nickel sheet pin from B + and B-with shear pliers, and take out the battery.

- Before the waste appliance, the battery must be removed from the appliance;
- When taking out the battery, the appliance must be cut *off*;
- The battery shall be disposed of safely.

This device contains a non-replaceable battery. The equipment and battery must be disposed of properly when the battery life is over or leaks.

● Size of product

Product name	Oclean X Pro 20 WiFi Smart Sonic Electric Toothbrush
Product model	NIB04
Rated input	5V=0.5A*
Rated input power	2.5W
Battery capacity	1100mAh
Compatible devices	Android with compatibility5.0,IOS9.0And above, with Bluetooth support4.2
Frequency range	2400-2483.5MHz
Waterproof grade	IPX7
Charging method	Type-C

*When charging, an adapter with an output of 5V=0.5A (the maximum power of the output is less than 15W) and in compliance with the requirements of IEC 60335-1 must be used for charging

Manufacturer: Shenzhen Yunding Information Technology Co., Ltd.

Manufacturer Address: 28G, Building 3, Dachong Business Center(Phase 3), No.18, Dachong 1st Road, Dachong Community, Yuehai Street, Nanshan District,Shenzhen City, Guangdong Province, China

EU Authorized Representative European Address: Krokuvos
str. 13, Agency Name: Mbrush UAB 09314 Vilnius, Lithuania

Email:
service@mobalt.eu Post
Code: LT44237