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
Oclean X10 Smart Sonic Electric Toothbrush User Manual

Safety Tips

This product is for 12 years old and above.
The product is compatible with Oclean replacement brush heads only. Oclean bears no responsibility for any consequences resulting from the misuse of non-certified third-party brush heads.
This product is for personal cleaning only and not for medical purposes or multi-user.
This product is for indoor use only. Do not place the toothbrush under direct sunlight for a prolonged time. It is recommended to be placed in dry well-ventilated areas with normal temperature.
Please keep the toothbrush handle dry or wipe it with a soft cloth before charging. After charging, please resume the waterproof charge port cover back to its original place and ensure that it is sealed. Oclean is not responsible for any consequences caused by improper operation.
This product can only be charged via the corresponding voltage demonstrated on the product.
The built-in battery can only be changed by trained professionals.

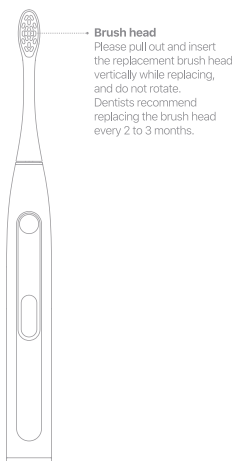Instruction

Up/down control
The toothbrush is operated via up-down control, tap the up button to switch modes, and the down button to switch intensities. There are 5 modes to select: Clean Mode, White Mode, Massage Mode, Sensitive Mode, and Polishing Mode. Each has 5 levels of intensities.

ON/OFF control
Once toothbrush is turned on, tap the ON/OFF button to start brushing your teeth. Tap the button to pause during the brushing session and then tap again within 15 seconds to resume its operation; Whilst in pause with over 15 seconds will inactivate the toothbrush. The toothbrush will be in hibitnate with a tap on the ON/OFF button. To reset the language setting, please press the "ON/OFF control" button for 8 seconds to revert to the factory setting.

Brushing program
White Mode is the default mode with 2 mins timer, and the toothbrush will remind the user to switch the brushing zone every 30 seconds by vibrating.

Smart brushing feedback

After each brushing session, there will be emoji feedbacks of the brushing report with a stars score rating. When score 3 stars, there will be a smiley emoji to demonstrate that you brush like a dentist. When scored less than 3 stars, there will be a tooth model to demonstrate areas that you have been missing while brushing. Then a supplementary brushing plan with according time will be generated for you.

Then simply tap the up/down control to activate the supplementary brushing plan. When completed, the toothbrush will show an emoji on the display, and end this brushing session.

Toothbrush charging

This product adopts the USB Type-C charging port. When charging is in need, remove the waterproof charge port cover and connect the charging cable to a power plug. When fully charged, please unplug the charging cable resume the cover. (as shown in the figure)

Overcharging may lead to premature battery depletion, thus affecting the user experience.

Thanks for using Oclean Smart Sonic Electric Toothbrush.

Note

Attention should be paid to those with the following conditions of oral sensitivities

Those who have recently undergone periodontal surgery or are undergoing Oral-related treatments.
Individuals with other symptoms orally include those with unstable dentures or have heavy gum bleeding.
Other individuals with sensitive and fragile oral issues.
This product is not recommended for those with the above conditions and should consult a professional dentist before use.
This product is not suitable for those with disabilities, mental illness, or lack of relevant experience and knowledge [including children (0-12 years old)] unless under the supervision of their guardians.

Low Battery Reminder

When the battery percentage is below 20%, a reminder will pop up on the display of your toothbrush from time to time to indicate its battery situation. When the battery is extremely low with less than 1%, your toothbrush and its screen will be off, and the motor will vibrate to alert you to charge the battery.

How to Charge

When charging is needed, remove the protection cover at the bottom of the toothbrush, connect the charging cable to the power plug and the toothbrush. During charging, the screen on the toothbrush will display a charging interface. Generally, it takes about 3 hours to fully charge the device.

Recommended tooth brushing method

Position the toothbrush at 45° between the bristles and gingiva line. It isn't necessary to press hard or scrub, just simply guide the brush


Start with the outside surfaces of the teeth. Guide the brush head slowly from tooth to tooth, holding the brush head in place for a few seconds

Mouthwash with water afterward.

Device cleaning and maintenance

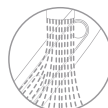Toothbrush handle and brush head cleaning
Direct rinse under running water with a cover on the brush bristles. After washing, air dries the brush head, and no need to pull the bristles.

Product Specifications

Product Name	Oclean X10 Smart Sonic Electric Toothbrush
Model	R3100
IP Rating	IPX7
Input Power	2.5W
Charging	Type-C charging
Input	5V= 0.5A
Battery life*	60 days
Battery Capacity	1200 mAh
Charging time	3 hours
Implementation standard	GB 4706.1-2005,GB 4706.59-2008

* Calculated under whitening mode, 2 minutes/time, 2 times/day.

This device complies with Part 15 of the FCC rules. Operation is subject to the following two conditions:
1) this device may not cause harmful interference, and
2) this device must accept any interference received, including interference that may cause undesired operation

EU Authorized Representative
Agency name: Mbrush UAB
Representative: Algirdas Tevelis
European Address: Totori g.3, Kaunas, Lithuania
Post Code: LT44237
Email: service@mobalt.eu
Phone: +37067052772

-
1. Do NOT disassemble, strike, extrude or throw the device into a fire.
 2. Do NOT use the toothbrush if its battery is bulging or has expanded.
 3. Do NOT place in a high-temperature environment.
 4. Do NOT use the device if its battery has been subjected to immersion in water.
 5. Before the appliance is discarded, its battery must be removed from it.
 6. Before the battery is removed, the appliance must be powered off.
 7. The battery shall be disposed of safely.

Name of parts	Toxic and harmful substance or element							Environment Protection
	Pb	Hg	Cd	Cr6+	PBBs	PBDEs	4 Phthalates	
PCBA	○	○	○	○	○	○	○	10 years
Plastic parts	○	○	○	○	○	○	○	
Metal hardware	○	○	○	○	○	○	○	

Name of parts	Toxic and harmful substance or element							Environment Protection
	Pb	Hg	Cd	Cr6+	PBBs	PBDEs	4 Phthalates	
Vibrating motor	○	○	○	○	○	○	○	10 years
Accessories	○	○	○	○	○	○	○	
Battery	○	○	○	---	---	---	---	

Remarks

- This indicates that the content of toxic and hazardous substances in all homogeneous materials of this component is less than the limited quantity as specified in RoHS 2.0 Directive: 2011/65/EU & (EU) 2015/863.
- × This indicates that the content of toxic and hazardous substances in at least one of the homogeneous materials of this component is more than the limited quantity as specified in RoHS 2.0 Directive: 2011/65/EU & (EU) 2015/863.
-
- * Except for very few materials, due to technical development reasons, most of the materials used are environmentally friendly. Please use this product with confidence during the environmental protection period.

Guarantee Certificate of Oclean Smart Sonic Electric Toothbrush

User Information

Name: _____ Phone: _____ E-mail: _____

Address: _____ ZIP code: _____

Retailer Information

Retailer: _____ Phone: _____

Address: _____ ZIP code: _____

Invoice no.: _____ Sales date: _____

Product Information

Product type: _____ Manufacturing site: _____

Warranty Card of Oclean Smart Sonic Electric Toothbrush

Delivery date	Fault description/reason/disposal	Equipment employed	Acceptance date	Serviceman

Please request a maintenance report from the maintenance center and keep it together with the warranty certificate. (See the page on the left.)

Shenzhen Yunding Information Technology Co., Ltd., in strict accordance with the Responsibility Regulations on Maintenance, Replacement and Sales Return of Household Appliances, offers a guarantee for consumers who purchase the Oclean Smart Sonic Electric Toothbrush. In the case of any inconsistency between the service articles and the provisions of guarantee, the provisions of the guarantee and the contents of the guarantee certificate shall govern.
If the product you have purchased is found to have a fault or defect, please contact our service line at 4008-365-909.

Manufacturer: Shenzhen Yunding Information Technology Co., Ltd.
Manufacturer's address: 28G, Building 3, Dachong Business Center (phase III), No.18 Dachong 1st Road, Dachong Community, Yuehai Street,Nanshan District, Shenzhen, Guangdong, China.
Production address: No.89, Dingshan Middle Road, Dongfushan Street, Haicang District, Xiamen, Fujian, China
Made in China

Website: www.oclean.com

IPX7 CE FCC IC EAC